


# Flexible Illumination Systems for Optimal Visualisation

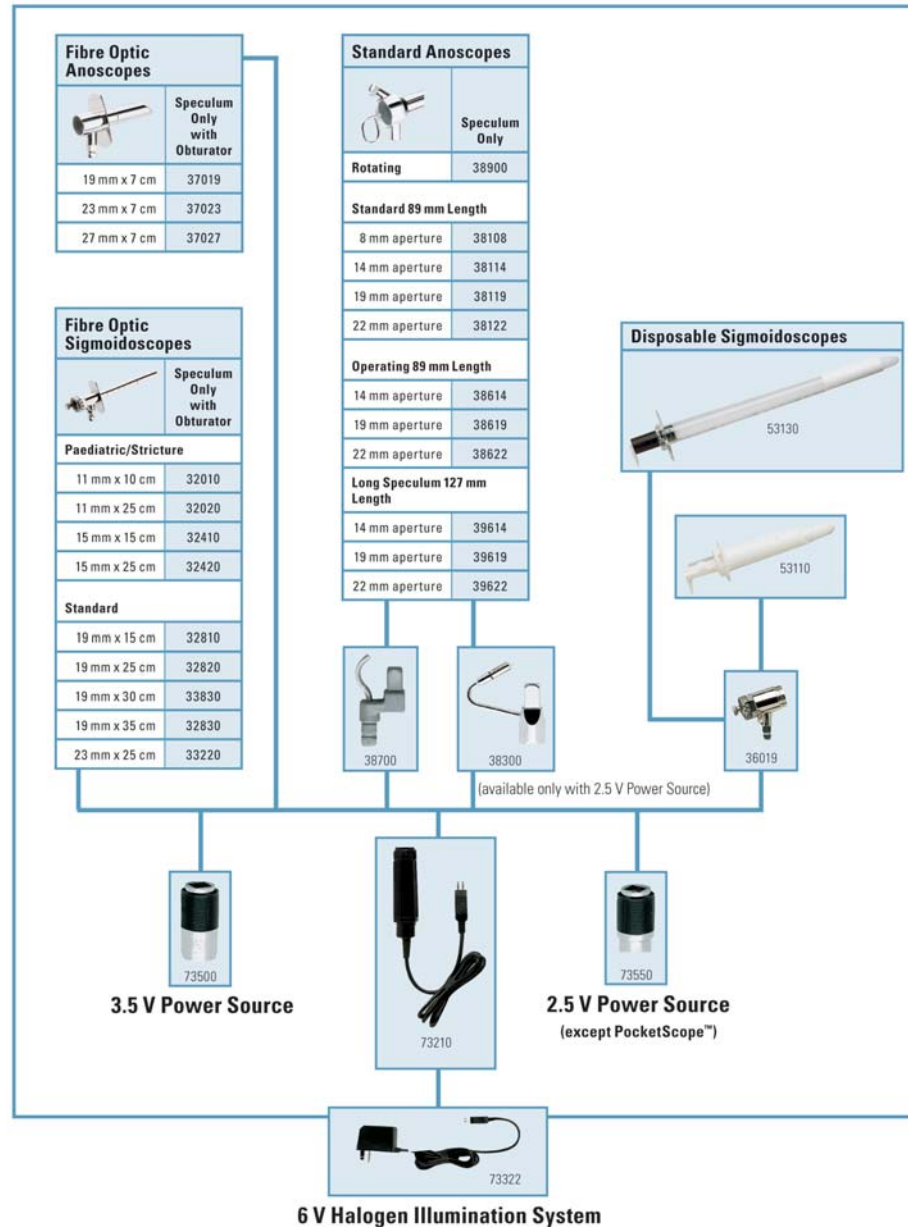

## Sigmoidoscope/Anoscope Lighting Systems

Powerful illumination systems.

- > Halogen light provides consistent, long-lasting illumination (73500, 73550 only)
- > Fibre optics provide a cool distal ring of light with no reflections or obstructions (36019, 38700 only)
- > Sigmoidoscopes/anoscopes can be used with a choice of illumination systems—a 6 V system with optional on/off switch and rheostat control, or with Welch Allyn's 3.5 V/2.8 V or 2.5 V hand-held power sources (with adaptors 73500 and 73550).
- > 73500 adapter fits most Welch Allyn 3.5 V power handles
- > 73550 adapter fits most Welch Allyn 2.5 V power handles (except PocketScope™)

## Sigmoidoscope/Anoscope Lighting Systems and Accessories

<b>36019</b>	Fibre Optic Light Head for Disposable Sigmoidoscopes and Anoscopes	<b>73550</b>	2.5 V Halogen Handle Adapter
<b>38300</b>	2.5 V Vacuum Light Carrier	<b>30200</b>	Insufflation Bulb, Complete
<b>38700</b>	Fibre Optic Light Carrier for use with Standard Anoscopes and Light Adapter	<b>30130</b>	30 cm Suction Tube (Standard)
<b>73210</b>	Light Handle with 0.9 m cord	<b>30140</b>	40 cm Suction Tube
<b>73308</b>	Optional On/Off Switch	<b>34000</b>	Rectal Probe
<b>73322</b>	Transformer with 1.5 m cord (EUR)	<b>34100</b>	Rectal Hook
<b>73400</b>	Optional Rheostat Control	<b>00200-U</b>	2.5 V Vacuum Lamp (small) (for 38300)
<b>73500</b>	3.5 V Halogen Handle Adapter	<b>03100-U</b>	3.5 V Halogen Lamp (for 73500)
		<b>03400-U</b>	2.5 V Halogen Lamp (for 73550)
		<b>07800-U</b>	6 V Halogen Lamp (for 73210)

For more information, contact us at +353 46 906 7790 today!

Welch Allyn European Customer Service, Navan Business Park,  
Dublin Road, Navan, County Meath, Ireland  
Tel: +353 46 906 7790 Fax: +353 46 906 7754  
www.welchallyn.com

Welch Allyn PTY. LTD, S. Africa  
394 Surrey Ave, Ferndale, Randburg 2194 South Africa  
P O Box 809, Fontainebleau, South Africa 2032  
Tel: 0 11 777 7555 Fax: 0 11 777 7556  
www.welchallyn.com

# Quality and Value in Rigid Endoscopy

Welch Allyn offers a wide range of rigid endoscopy products for anoscopy and sigmoidoscopy, complemented by versatile fibre optic illumination systems for optimum visualisation.

The sigmoidoscope range includes both reusable and disposable instruments and a variety of sizes to accommodate differing requirements. Welch Allyn fibre optic anoscopes represent the finest in design and craftsmanship for optimal performance.



## Welch Allyn Disposable Sigmoidoscopes

One of our most popular designs, the Welch Allyn Disposable Sigmoidoscope may be used in conjunction with either a 6 V illumination system or traditional Welch Allyn 3.5 V power handles.

- > Highly visible graduations for greater accuracy and control
- > Glare-free frosted coating ensures maximum light output at distal tip
- > Water-activated lubricating gel to facilitate easy insertion
- > Individually wrapped for safety and convenience
- > Fibre-optic illumination system provides a cool 360 degree ring of light
- > Hinged window with neoprene seal to avoid fogging during insufflation

<b>53130</b>	KleenSpec® Disposable Sigmoidoscope with Obturator (case of 100), (add 'L' suffix for lubricated version)
<b>73500</b>	3.5 V Halogen Handle Adapter
<b>71020-A</b>	3.5 V Rechargeable Handle
<b>73210</b>	Light Handle with 0.9 m Cord
<b>36019</b>	Fibre Optic Light Head
<b>73322</b>	Transformer with 1.5 m Cord (EUR)

## Insufflation Systems

An insufflation system can be used with Welch Allyn sigmoidoscopes to ease insertion of the instrument. Hydrophobic filters are available to prevent contamination of the insufflation bulb. They are attached using a Luer lock fitting.

<b>30200</b>	Insufflation Bulb, complete
<b>30215</b>	Luer Lock for use with hydrophobic filters (pack of 10)
<b>31210</b>	Hydrophobic filters (pack of 10)



## Standard and Paediatric/Stricture Sigmoidoscopes

All Welch Allyn reusable sigmoidoscopes feature:

- > Autoclavable stainless steel for convenience and durability
- > Hinged window with a neoprene seal to avoid fogging during insufflation
- > Graduated tube with removable tip for fibre cleaning (metal only)



### Standard Sigmoidoscopes

<b>32810</b>	19 mm x 15 cm Standard Sigmoidoscope, Speculum only with Obturator
<b>32820</b>	19 mm x 25 cm Standard Sigmoidoscope, Speculum only with Obturator
<b>33830</b>	19 mm x 30 cm Standard Sigmoidoscope, Speculum only with Obturator
<b>32830</b>	19 mm x 35 cm Standard Sigmoidoscope, Speculum only with Obturator
<b>33220</b>	23 mm x 25 cm Standard Sigmoidoscope, Speculum only with Obturator



### Fibre Optic Sigmoidoscopes

<b>32010</b>	11 mm x 10 cm Fibre Optic/Stricture Speculum only with Obturator
<b>32020</b>	11 mm x 25 cm Fibre Optic/Stricture Speculum only with Obturator
<b>32410</b>	15 mm x 15 cm Fibre Optic/Paediatric Speculum only with Obturator
<b>32420</b>	15 mm x 25 cm Paediatric/Stricture Speculum only with Obturator

## Welch Allyn Anoscopes

A complete line of high-performance and reusable anoscopes for safe and efficient examinations.



### Disposable Anoscope

<b>53110</b>	KleenSpec® Disposable Anoscope with Obturator (box of 25)
--------------	---

### Disposable Anoscope Accessories

<b>73210</b>	Light Handle with 0.9 m Cord
<b>73222</b>	Light Handle (73210) plus Transformer (73322)
<b>36019</b>	Fibre Optic Light Head
<b>73322</b>	Transformer with 1.5 m Cord (EUR)



### Rotating Anoscope

Anoscope Only 68 mm Speculum Length

<b>38900</b>	Rotating Speculum Anoscope with Obturator, 20 mm inside diameter, 22 mm outside diameter
--------------	--



### Fibre Optic Anoscope

Speculum Only with Obturator

<b>37019</b>	19 mm x 7 cm Fibre Optic Anoscope, Speculum Only with Obturator
<b>37023</b>	23 mm x 7 cm Fibre Optic Anoscope, Speculum Only with Obturator
<b>37027</b>	27 mm x 7 cm Fibre Optic Anoscope, Speculum Only with Obturator

*Welch Allyn reusable anoscopes are also available in Long-Speculum, Operating and Standard configurations. Please contact us for further information on sizes and product codes.*